
HeartSine®

Pad-Pak™ and Pediatric-Pak™

Set-up

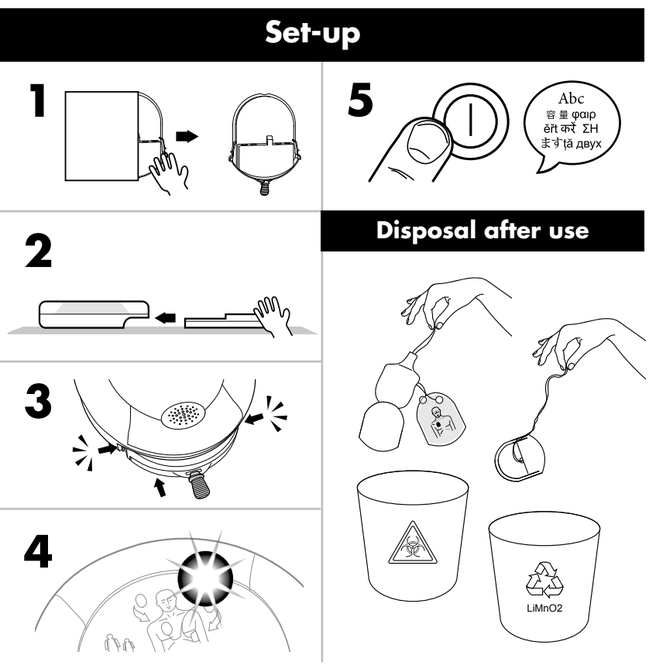Disposal after use

HeartSine

CE 0123

HeartSine Technologies, Ltd.
207 Airport Road West
Belfast
Northern Ireland
BT3 9ED
United Kingdom
Tel +44 28 9093 9400
Fax +44 28 9093 9401

Distributed in U.S. by:

Stryker Emergency Care
11811 Willows Road NE
Redmond, WA, 98052 U.S.A.
Toll Free 866 478 7463
Tel +1 215 860 8100
Fax +1 215 860 8192

heartsine.com
heartsinesupport@stryker.com


Stryker Australia Pty Ltd
8 Herbert St
St Leonards NSW 2065
Australia

Stryker European
Operations B.V.
Herikerbergweg 110
1101 CM Amsterdam
Netherlands

HeartSine samaritan PAD: UL Classified. See complete marking on product.

© 2022 HeartSine Technologies. All rights reserved.

Made in U.K.
Date of issue: 11/2022
H017-013-363-AF

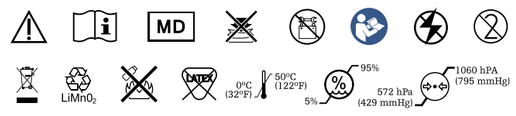EN-US (US only)

Set-up

1. Remove from packaging.
2. Place HeartSine samaritan PAD (AED) face up on a flat surface and slide in Pad-Pak or Pediatric-Pak.
3. Listen for "double click" to indicate it is properly inserted.
4. Verify that green status indicator is blinking, indicating the AED is ready for use.
5. Turn on AED. You should hear expected language and should not hear a warning message. Turn off AED.

Disposal after use

- Remove Pad-Pak. Cut wires. Stick pads together. Dispose of pads as infectious waste. Recycle lithium battery separately per local regulations

Pad-Pak and Pediatric-Pak contain:

- Pair of disposable single use defibrillation pads
- $LiMnO_2$ (18V – 1500mAh) non-rechargeable battery

Indications for use

The Pad-Pak and Pediatric-Pak, when used with a HeartSine samaritan PAD, are indicated for use on victims of cardiac arrest. See the AED user manual or user handbook for the exhibiting signs of cardiac arrest specific to each device.

- **Pad-Pak:** Intended for use on adults or children > 8 years or > 25 kg (55 lb)
- **Pediatric-Pak:** Intended for use on children between the age of 1-8 years or < 25 kg (55 lb)

Battery capacity

>60 shocks or 6 hours battery use.

Warnings

DO NOT USE if expiry date is exceeded.

DO NOT USE if the pouch is opened or damaged. This may result in the electrode gel being dry.

Cautions

Single use ONLY. Reuse may cause AED to be unable to deliver therapy leading to a failure to resuscitate. It may also lead to cross-infection between patients.

DO NOT PULL the green tab on the Pad-Pak during set-up.

If you hear a warning message when you turn on your HeartSine samaritan PAD, consult the troubleshooting section of the HeartSine samaritan PAD user manual or user handbook.

For full instructions for use, indications for use, contraindications and safety information, and disposal instructions, refer to your HeartSine samaritan PAD user manual or user handbook. Each can be found at heartsine.com.

Notes

- Pediatric-Pak is only suitable for use with HeartSine samaritan PADs with the symbol
- The TSO/ETSO-certified Pad-Pak must be used on board commercial fixed wing aircraft
- DO NOT STORE the Pediatric-Pak next to magnetically sensitive storage media as it contains a magnetic component (surface strength 6500 gauss)

Please report any serious incident that occurs with this device to HeartSine Technologies, Ltd. and to the U.S. Food and Drug Administration as per local regulations.

EN-UK

Set-up

1. Remove from packaging.
2. Place HeartSine samaritan PAD (AED) face up on a flat surface and slide in Pad-Pak or Pediatric-Pak.
3. Listen for "double click" to indicate it is properly inserted.
4. Verify that green status indicator is blinking, indicating the AED is ready for use.
5. Turn on AED. You should hear expected language and should not hear a warning message. Turn off AED.

Disposal after use

- Remove Pad-Pak. Cut wires. Stick pads together. Dispose of pads as infectious waste. Recycle lithium battery separately per local regulations

Pad-Pak and Pediatric-Pak contain:

- Pair of disposable single use defibrillation pads
- $LiMnO_2$ (18V – 1500mAh) non-rechargeable battery

Intended purpose

The Pad-Pak or Pediatric-Pak accessory contains the battery to power the HeartSine samaritan PAD and two electrode pads to provide the electrical connection to the patient's chest for delivery of defibrillation shock. The delivery of a defibrillation shock can stop disruption to the heart's normal rhythm and restore blood flow.

Indications for use

The Pad-Pak and Pediatric-Pak, when used with the HeartSine samaritan PAD, are indicated for use on victims of cardiac arrest who are exhibiting the following signs: unconscious, not breathing and without circulation.

- **Pad-Pak:** Intended for use on adults or children > 8 years or > 25 kg (55 lb)
- **Pediatric-Pak:** Intended for use on children between the age of 1-8 years or < 25 kg (55 lb)

Contraindications

Do not use if patient is responsive or conscious.

Battery capacity

>60 shocks or 6 hours battery use.

Warnings

DO NOT USE if expiry date is exceeded.

DO NOT USE if the pouch is opened or damaged. This may result in the electrode gel being dry.

Cautions

Single use ONLY. Reuse may cause AED to be unable to deliver therapy leading to a failure to resuscitate. It may also lead to cross-infection between patients.

DO NOT PULL the green tab on the Pad-Pak during set-up.

If you hear a warning message when you turn on your HeartSine samaritan PAD, consult the troubleshooting section of the HeartSine samaritan PAD user manual or user handbook.

For full instructions for use and safety information refer to your HeartSine samaritan PAD user manual or user handbook. Each can be found at heartsine.com.

Notes

- Pediatric-Pak is only suitable for use with HeartSine samaritan PADs with the symbol
- The TSO/ETSO-certified Pad-Pak must be used on board commercial fixed wing aircraft
- DO NOT STORE the Pediatric-Pak next to magnetically sensitive storage media as it contains a magnetic component (surface strength 6500 gauss)

Please report any serious incident that occurs with this device to HeartSine Technologies, Ltd. and to your national competent authority or other local regulatory authority as per local regulations.

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